


NEGATIVE PRESSURE PUMP

# NewEraII®

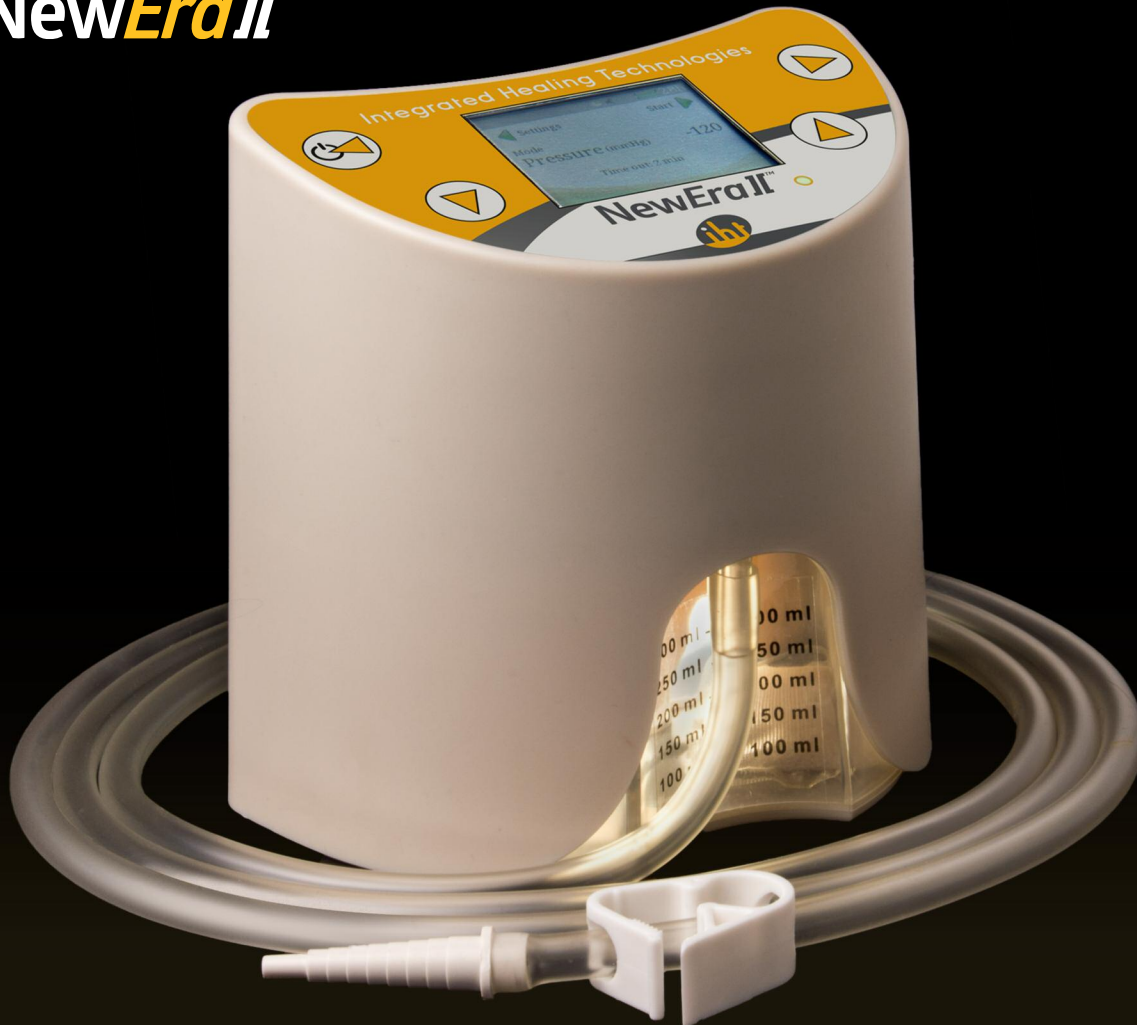

## NewEraII® CANISTER

MODEL NUMBER: 16690-NEWERA

- Light weight
- Ergonomic design
- Fits comfortably on hip
- Carrying bag included
- Pre-attached tubing and nozzle with tubing clamp
- Inconspicuously located canister
- Single touch canister removal
- Reliable pump motor





The light weight of the NewEra-II with its ergonomic design makes this pump ideal for the ambulatory setting. The pump is supplied with a carrying bag and has a concave shape to fit comfortably on the hip. The canister is inconspicuously and conveniently located within the exterior housing of the pump. The motor that operates the pump to deliver negative pressure and remove fluids is the same reliable motor that exists in IHT's NewEra-I NPWT Pump.



### NewEra II PUMP SPECIFICATIONS

| REQUIREMENT                  | SPECIFICATION            | NOTES  |
|------------------------------|--------------------------|--|
| Portable                     | Yes                      | Carry bag available  |
| Weight                       | 1.82 lbs.                | With empty canister  |
|                              | 2.48 lbs.                | Maximum weight with full canister  |
| Pressure Range               | -10 to -200 mmHg         | Adjustable in 5 digit increments   |
| Battery Life                 |                          | Up to 10 hours   |
| Cycles                       | Yes                      | Intermittent or Continuous   |
| Displays and alarm activated | Digital                  | Displays pressure setting, mode of operation                                 |
|                              | Status LED               | Operational LED  |
|                              | Other LED                | Slow Mode controls rate of negative pressure increase or decrease            |
| Auto Lockout                 | Yes                      | User controlled time settings  |
| Blue Tooth Data Capture      | Yes                      | Alarms, On/Off, runtimes, notifications                                      |
| Pressure Tolerance Alarm     | Yes                      | If the vacuum is from 5 to 80 mmHg below the user controlled set point       |
| High Vacuum Alarm            | Yes                      | If the vacuum is from 5 to 80 mmHg higher than the user controlled set point |
| Other Alarms canister full   | Yes                      | Low battery, Tube blockage, Inactive pump,                                   |
| Alarm Mute Button            | Yes                      |  |
| Main Power Option            | Yes                      |  |
| Power Input                  | Internal battery or 100V |  |

### NewEra II CANISTER SPECIFICATIONS

| REQUIREMENT                   | SPECIFICATION | NOTES  |
|-------------------------------|---------------|--|
| Canister Sizes                | 300ml         |  |
| Solidifier with Canister      | Yes           |  |
| Fluid Measurement Marks       | Yes           |  |
| Single touch canister removal | Yes           |  |
| Antibacterial Filter          | Yes           | 3 micron bacterial filter with a 32 mm diameter to prevent microbial contamination |
| Tubing Clamps                 | Yes           |  |

**Integrated Healing Technologies™**